



**BD Luer-Lok™  
Tip**

Generally used for injections requiring a secure connection of the syringe to another device. The tip is threaded for a 'locking' fit, and is compatible with a variety of needles, catheters, and other devices.

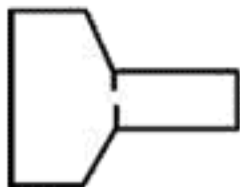


**Luer Slip Tip**

A friction-fit connection that requires the clinician to insert the tip of the syringe into the needle hub or other attaching device in a push-and-twist manner. This will ensure a connection that is less likely to detach. Simply sliding the attaching device onto the syringe tip will not ensure a secure fitting.

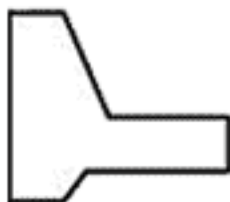

**Eccentric Luer  
Slip Tip**

Allows for work requiring closer proximity to the skin. Generally used for venipunctures and aspiration of fluids. (Also see luer slip instructions above).

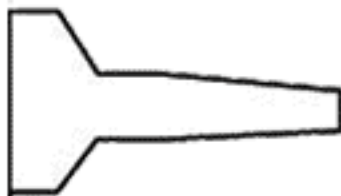

**Catheter Tip**

Used for flushing (cleaning) catheters, gastrostomy tubes and other devices. Insert catheter tip securely into catheter or gastrostomy tube. If leakage occurs, refer to your facility's guidelines.