


**BD Luer-Lok™
Tip**

Generally used for injections requiring a secure connection of the syringe to another device. The tip is threaded for a 'locking' fit, and is compatible with a variety of needles, catheters, and other devices.

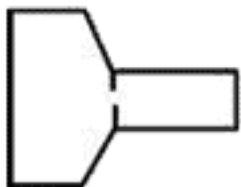


Luer Slip Tip

A friction-fit connection that requires the clinician to insert the tip of the syringe into the needle hub or other attaching device in a push-and-twist manner. This will ensure a connection that is less likely to detach. Simply sliding the attaching device onto the syringe tip will not ensure a secure fitting.



**Eccentric Luer
Slip Tip**

Allows for work requiring closer proximity to the skin. Generally used for venipunctures and aspiration of fluids. (Also see luer slip instructions above).

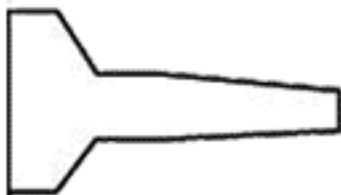

Catheter Tip

Used for flushing (cleaning) catheters, gastrostomy tubes and other devices. Insert catheter tip securely into catheter or gastrostomy tube. If leakage occurs, refer to your facility's guidelines.